



Device and humidifier overview

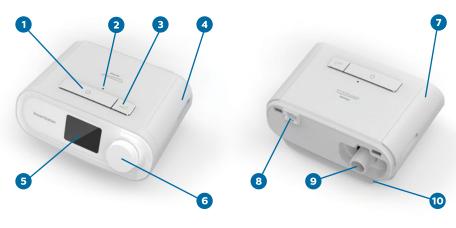

Device feature

_			
1.	Therapy On/Off button	6.	Control dial
2.	Ambient light sensor	7.	Door, accessory access
3.	Ramp button	8.	Humidifier connector
4.	Door, SD card & filter access	9.	Air outlet port
5	Display screen	10	Power inlet

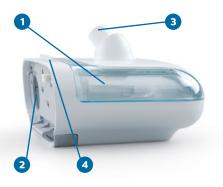



Device feature

1.	Humidifier	4.	Humidifier lid release latch
2.	Air inlet port	5.	Water tank
3	Air outlet port	6.	Maximum fill lines

Setup

Your DreamStation sleep therapy system comes with the device, power supply with AC cord, flexible tube and an optional heated humidifier.



Apply power

Plug one end of the AC cord into the power supply, and the other end into an electrical outlet. Plug the green power supply connector into the green power inlet on the therapy device.



Access the humidifier water tank

Slide the humidifier lid release latch toward the therapy device and lift up on the humidifier lid until it is completely open.



Fill the humidifier water tank

Lift the tank out of the humidifier base. Rinse the tank with tap water before your first use. Then fill the tank with distilled water up to the maximum fill line.



Replace water tank

Carefully replace the water tank into the humidifier and close the humidifier lid until it snaps shut.



Attach the tube and mask

Attach the flexible tubing to the outlet port on the top of your humidifier, or on the back of your therapy device if you are not using a humidifier.

Connect the other end of the tubing to your mask and headgear.



Starting/stopping therapy

To start or stop therapy, press the therapy on/off button.

Caution: US federal law restricts these devices to sale by or on the order of a physician. Warning: see user manual for complete instructions.

Navigating in the menu

When therapy is on:



Therapy pressure screen

Device feature

- Therapy pressure
- 2. Adjustable tube temperature setting
- 3. Adjustable humidifier setting
- 4. Enabled features

Humidifier settings

While the device is delivering therapy, you can adjust tube temperature or humidifier settings. Rotate the control dial to choose either setting. Press and rotate the dial to change the setting.

When therapy is off: in the patient menu, you can choose between the following screens. Only the menus available and enabled on your device will display.





My Info

This menu provides summary statistics of your therapy use.



Preheat (if available)

This function lets you warm up your humidifier for 30 minutes before starting a therapy session.



My Provider

This menu contains information that your provider may direct you to read to them so they can better assist you over the phone.



My Setup

This menu contains comfort settings that you can adjust as needed.

Caring for your DreamStation

Cleaning your DreamStation

DO NOT immerse the device in any fluids.

- Unplug the device, then wipe the outside of the device with a cloth slightly dampened with water and a mild detergent.
- Let the device dry completely before plugging in the power cord.

Cleaning the tubing

Clean the flexible tubing before first use and daily.

- · Disconnect the flexible tubing from the device.
- Gently wash the tubing in a solution of warm water and a mild detergent. Rinse thoroughly. Air dry. Inspect the tubing for damage or wear.
- Refer to the humidifier manual for the instructions on how to clean the optional heated tube.

Cleaning or replacing the filters

Disconnect the device from the power source.

- · Remove the filter(s) from the device.
- To clean the reusable blue pollen filter, detach the light-blue ultra-fine filter if applicable, and set aside. Take the reusable filter to the sink, turn it upside down, and run warm tap water through the white filter media to rinse away any debris. Allow the filter to air dry completely before reinstalling.
- DO NOT clean the ultra-fine filter. If the light-blue ultra-fine filter is dirty or torn, replace it.
- · Reinstall the filters.





Mask options

A comfortable, well fitted mask is one of the key factors to therapy success.

With innovative cushion designs for excellent seal and comfort, our Philips Respironics masks help CPAP users get the sleep they deserve.



Nasal masks

Nasal masks are generally smaller and lighter than full face masks and may be better for higher air pressure than minimal contact masks.



Pillows masks

These masks may be a good solution if you feel claustrophobic, if you are an active sleeper, or if you have facial hair that interferes with other masks.



Full face masks

This mask may be ideal if you suffer from nasal obstruction or frequent congestion from allergies or cold symptoms. It's also suitable if you continue to breathe through your mouth at night.



DreamWear Nasal mask



DreamWearGel pillows mask



DreamWearFull face mask*



Minimal contact nasal mask



Nuance and Nuance Pro Pillows mask



Amara View Full face mask

For more information contact your provider or visit **sleepapnea.com**

DreamStation accessories



PAP battery kit

Lithium ion battery with integrated uninterruptible power supply (UPS) for DreamStation therapy device.



DreamStation travel case

The DreamStation travel case* gives patients a convenient, compact, and discreet way to travel with their DreamStation sleep therapy system and supplies.



PAP travel briefcase

The PAP travel briefcase** gives patients an all-in one carrying option that accommodates most PAP therapy systems and humidifiers, a laptop, and other small items.

For more information contact your provider or visit sleepapnea.com/products/accessories/travel-accessories

DreamStation Goportable travel CPAP machine



DreamStation Go and small travel kit

DreamStation Go is an indispensable travel partner, packing powerful performance into a light, streamlined package. It even has a convenient USB port to charge your mobile device. Also available is a small travel (as shown) or medium travel kit.

Optional overnight battery

Compact and lightweight with a recharge time of just 4-5 hours.



For more information contact your provider or visit sleepapnea.com/dreamstationgo

^{*}DreamStation, mask, and tubing not included. The case accommodates the Philips Respironics DreamStation sleep therapy system (blower only, no humidifier), the patient's mask and tubing, a 65W power supply w/ USB port, and a removable luggage strap.

^{**}DreamStation, power supply and tubing not included. The briefcase accommodates most PAP therapy devices and accessories and standard-size laptop computers.

Your

replacement schedule

The replacement schedule provided below are general guidelines for PAP therapy equipment. Everyone's insurance replacement schedule will vary. Please check with your home care provider to determine what your insurance coverage allows.

Equipment	Schedule	My insurance plan
Mask	1 per 3 months	
Nasal replacement cushion	2 per month	
Pillows replacement cushion	2 per month	
Full-face cushion	1 per month	
Headgear	1 per 6 months	
Tubing (standard or heated)	1 per 3 months	
Filter disposable	2 per month	
Filter reusable	1 per 6 months	
Water chamber	1 per 6 months	
Therapy device	1 per 5 years	

Keep a recordof your PAP information

Your provider name:	_Number:				
Device name and serial number:					
Humidifier serial number:					
PAP settings:Flex (if applicable):	Ramp:				
Mask name and size:					
Setup date:					

© 2018 Koninklijke Philips N.V. All rights reserved. Specifications are subject to change without notice. Trademarks are the property of Koninklijke Philips N.V. or their respective owners.

www.philips.com/respironics

Respironics Inc. 1010 Murry Ridge Lane Murrysville, PA 15668 USA



1130615

1130615 ROO JDW XX/XX/2018